



Dimensions

(Foot Print): W x L x H

Foot Print 10" x 22" x 8"
25.4 cm x 55.9 cm x 20.3 cm

Weight: 50 lb
22.7 Kg

Electrical Components Specifications

Motor RPM max. 12,000 RPM
an Power 208-230 V, 1 Phase, 20 A
Requirements:

max. 6,000 RPM
115V, 1 phase, 20 A

Motor Classification: DC Servo ^(a)
Brushless

Control Panel: Intrinsically Safe ^(a)

Materials of Construction Options (for wetted surfaces)

316L SS Hastelloy® C
Titanium Client's choice

Control Panel

The Magellan® system comes with an integrated electronic control panel and customizable software, data logging capability, and several external interfaces including RS-232, RS-422/485, USB and 10Base-T/100Base-TX Ethernet Port for LAN connection. A CompactFlash® socket is provided so that Flash cards can be used to collect trending and data logging information as well as to store larger configuration files.

Contact Kreido Laboratories regarding Heating & Cooling Requirements and support equipment. Working volumes are customized to each specific application.

Magellan®

Benchtop Quick Screening – External Heat Exchanger



The STT® system is a combination of patented processes and equipment for the manufacture of chemicals and pharmaceuticals. The Magellan™ Series of reactors is well suited for research and development labs.

The Magellan® type reactor is designed for early stage research where there is a need to quickly identify and confirm the validity of a reaction. The Magellan® series was developed to use small amounts of chemicals, so it is a good choice for applications where the reactants are expensive, such as pharmaceutical applications, or are hazardous (explosive or toxic). Experiments can be performed in a minimum amount of time and usually with only the simple support equipment typically found in most laboratories. Even though the reactor is small, it has been constructed using the same design principles as the Innovator® Series of pilot production reactors. Discoveries can be easily scaled to the larger equipment to supply larger amounts of product for testing or client evaluation.

On-line data logging is included in the system and can be connected to a company's internal customized software logging system. Customized connections can be provided to interface with real time analysis equipment such as HPLC or IR.

The STT® process and equipment are covered by several patents and patents pending worldwide. STT, Magellan, and Innovator are registered trademarks of Kreido Laboratories, all rights reserved. Kreido Laboratories is a trademark of Kreido Laboratories, all rights reserved. Hastelloy is a registered trademark of Haynes International, Inc. CompactFlash is a registered trademark of SanDisk Corporation. Copyright 2007 by Kreido.


Kreido Laboratories

1070 Flynn Road
Camarillo, CA 93021, USA


Office: 805-389-3499
Fax: 805-384-0989

Email: info@kreido.com
www.kreido.com

Magellan®

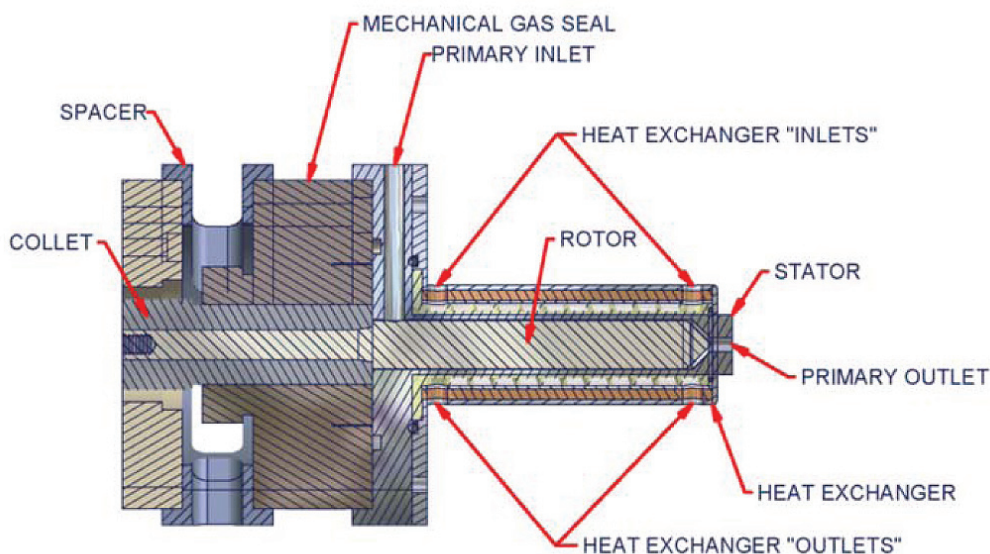

4 inch
reactor
length

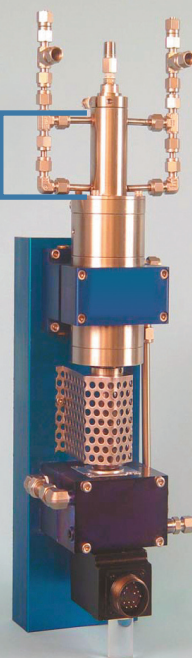

Design Specifications

Stator Inner Diameter (ID)	0.625", 15.9 mm		
Rotor Working Length	3.96", 10.1 cm		
Typical Rotor Outer Diameters (OD)	0.595" 15.1 mm	0.600" 15.2 mm	0.605" 15.4 mm
Gap	0.0150" 0.38 mm	0.0125" 0.32 mm	0.0100" 0.25 mm
Resulting Total Cavity Volume	2.06 ml	1.75 ml	1.43ml
Resulting Working Volume ^(e)	1.85 ml	1.54 ml	1.22 ml
Entry Port Configuration	2 ports standard, positioned at 180° with respect to each other, located in front of seal face ^(c)		
Exit Port Placement	Standard, at tip of rotor, as illustrated		

Operation Specifications

Rotation Rate	12,000 RPM maximum
Heat Exchanger – Required Pumping Volume	4 – 5.25 GPM (15-20 LPM) @ 5 PSIG

Example Shear Rates (sec⁻¹) for Stator ID = 0.625" at selected Rotor OD & RPM

Rotor OD	250 RPM	6,000 RPM	12,000 RPM
0.595" (15.1 mm)	519	12,462	24,923
0.600" (15.2 mm)	628	15,080	30,159
0.605" (15.4 mm)	792	19,007	38,013

Typical Seal Parameters

(continuous maximums)

Standard Seal

Gas Mechanical
300 PSIG
260°C
12,000 RPM

Optional Seals

Liquid Mechanical
Consult Factory

Note: Higher maximum, temperature, pressure, and rotor RPM values are interdependent.

The specifications provided here are general and Kreido Laboratories should be contacted for application recommendations when desired values fall out of the recommended range.

NOTES:

- Explosion Proof construction (with nitrogen purge) can be provided.
- Working Volume is based on the Rotor Working Length, the Rotor OD and the Gap where Gap = 1/2(Stator ID – Rotor OD). Custom Gap configurations can be provided but not with Gaps less than 0.010" or 0.25 mm
- More than two entry ports can be incorporated into the design (consult factory)

Kreido Laboratories

1070 Flynn Road
Camarillo, CA 93021, USA

Office: 805-389-3499
Fax: 805-384-0989

Email: info@kreido.com
www.kreido.com